


Fortress™

Reinforced Introducer Sheath System

Introducer sheath system with excellent kink resistance, easy insertion and superb radiopacity.

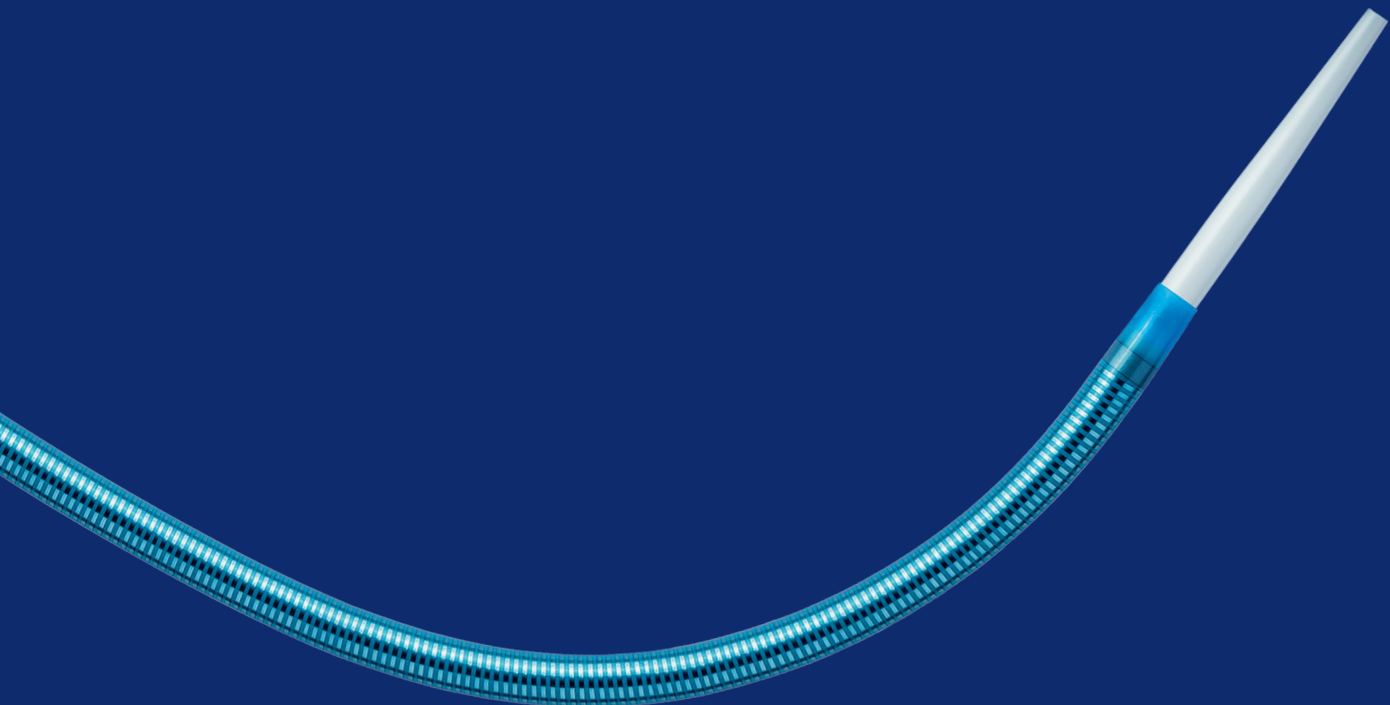

Fortress™ Introducer Sheath

Introducer sheath system with excellent kink resistance, easy insertion and superb radiopacity.


Excellent kink and deformation resistance


Coil reinforced sheath


The polymer embedded stainless steel coil shaft provides flexibility and excellent kink and deformation resistance.

PTFE liner

The non-ovaling characteristic combined with its friction lowering PTFE liner enables devices to easily pass through tight anatomical curves.

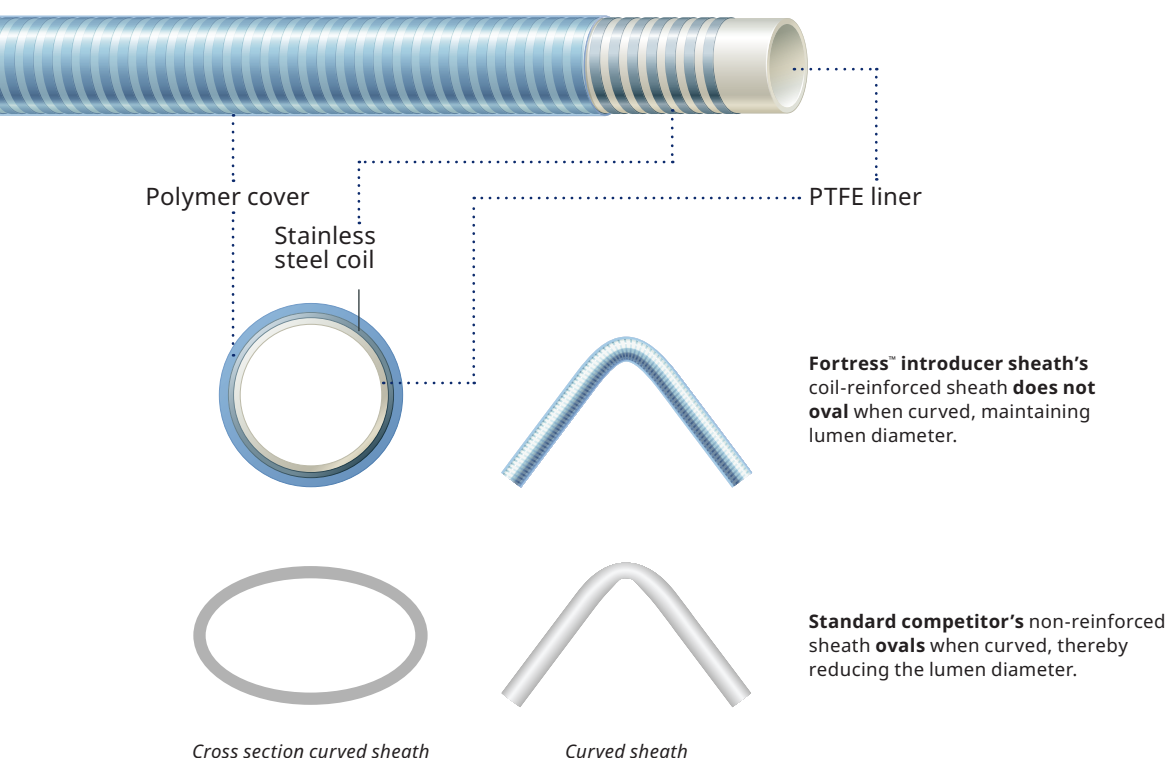

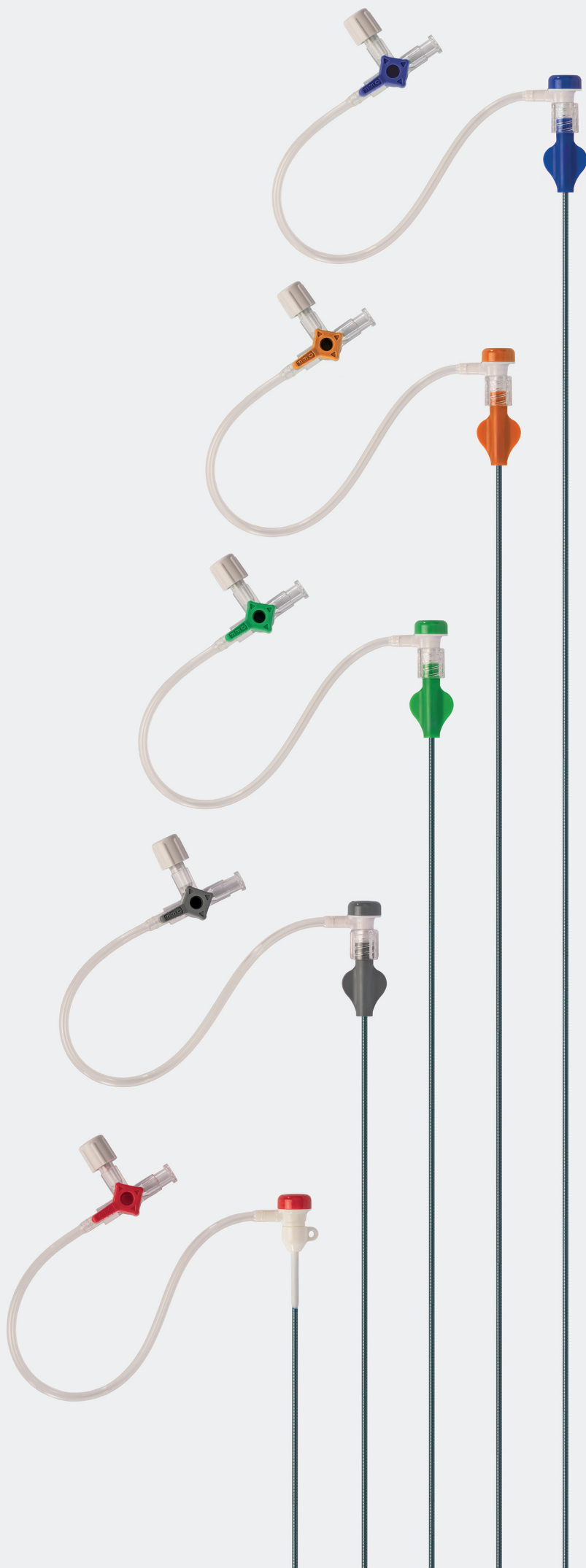

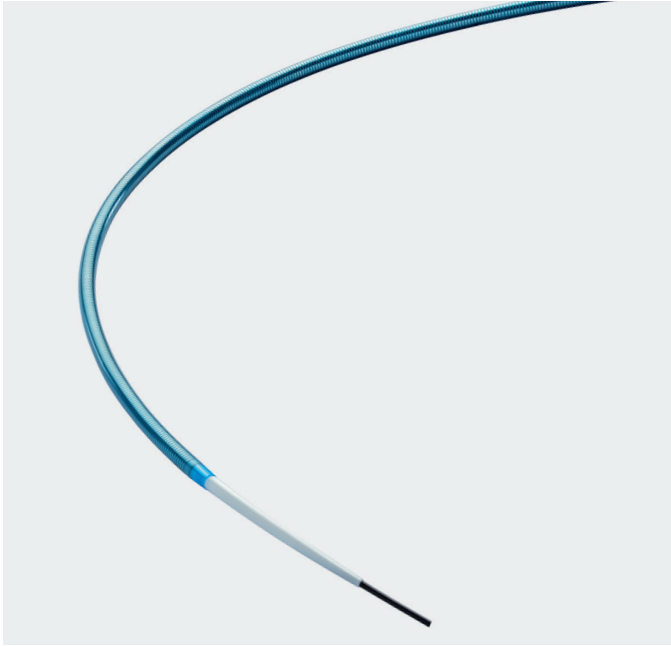

Easy insertion

Smooth tapered sheath and dilator tip

The sheath and dilator have a smooth tapered tip and an atraumatic tip design for ease of insertion.


Hydrophobic coated sheath¹

Best in class insertion performance² due to the hydrophobic sheath coating in combination with the smooth tapered sheath and dilator tip.

Superb radiopacity

Radiopaque sheath marker

The platinum alloy marker band provides superb tip radiopacity under fluoroscopic imaging.²

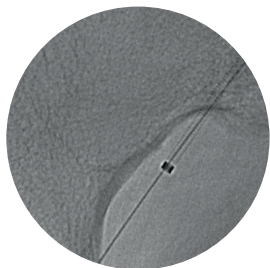


*Fortress™ 4F introducer sheath
in the pelvic area*

Minimal is optimal

4F

The 4F Solutions product portfolio allows for treatment of complex lower limb disease while reducing access site complications.

Fortress™ Introducer Sheath

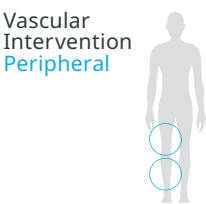

Indicated to provide access into the femoral, popliteal and infrapopliteal arteries.*

Technical data

INTRODUCER SHEATH	
French size	4F, 5F, 6F, 7F and 8F
Length	45, 65 and 90 cm
Sheath material	Stainless steel coil reinforced polymer, PTFE liner
Sheath tip	Smooth taper with radiopaque marker band located 3-5 mm from the sheath tip
Sheath shape	Pre-curved 45 cm Straight 45, 65 and 90 cm
Sheath coating	Hydrophobic, distal 30 cm (5F, 6F, 7F and 8F only)
DILATOR	
Recommended guide wire	0.035"
Tip	Smooth taper with radiopaque filler material
Shape	Pre-curved and straight
Hub	4 F: fixed hemostatic valve 5 F, 6 F, 7 F, 8 F: removable hemostatic valve Side-arm with color coded 3-way stop cock

Ordering Information

	SHEATH LENGTH	SHEATH SHAPE	INTRODUCER DESCRIPTION	CATALOGUE NUMBER	DILATOR(S) INCLUDED	UNITS PER BOX
4F	45 cm	Curved	Crossover	358813	Curved/Straight	5
	45 cm	Straight	Straight-45	358814	Curved/Straight	5
	65 cm	Straight	Straight-65	444486	Straight	5
	90 cm	Straight	Straight-90	444485	Straight	5
5F	45 cm	Curved	Crossover	386590	Straight	5
	45 cm	Straight	Straight-45	386591	Straight	5
	65 cm	Straight	Straight-65	444484	Straight	5
	90 cm	Straight	Straight-90	444483	Straight	5
6F	45 cm	Curved	Crossover	386593	Straight	5
	45 cm	Straight	Straight-45	386594	Straight	5
	65 cm	Straight	Straight-65	444482	Straight	5
	90 cm	Straight	Straight-90	444481	Straight	5
7F	45 cm	Curved	Crossover	452386	Straight	5
	45 cm	Straight	Straight-45	452382	Straight	5
	65 cm	Straight	Straight-65	452383	Straight	5
	90 cm	Straight	Straight-90	452388	Straight	5
8F	45 cm	Curved	Crossover	452387	Straight	5
	45 cm	Straight	Straight-45	452385	Straight	5
	65 cm	Straight	Straight-65	452384	Straight	5
	90 cm	Straight	Straight-90	452389	Straight	5

*Indication as per IFU.

References

- 1 5 F, 6 F, 7 F and 8 F.
- 2 Data on file.

Refer to IFU

Refer to QR code for a copy of the Instructions for Use for a complete listing of the indications, contraindications, warnings and precautions.



About Teleflex

As a global provider of medical technologies, Teleflex is driven by our purpose to improve the health and quality of people's lives. Through our vision to become the most trusted partner in healthcare, we offer a diverse portfolio with solutions in the therapy areas of anesthesia, emergency medicine, interventional cardiology and radiology, surgical, vascular access, and urology. We believe that the potential of great people, purpose driven innovation, and world-class products can shape the future direction of healthcare.

Teleflex is the home of Arrow™, Barrigel™, Deknatel™, LMA™, Pilling™, QuikClot™, Rüsch™, UroLift™ and Weck™ – trusted brands united by a common sense of purpose.

At Teleflex, we are empowering the future of healthcare. For more information, please visit teleflex.com.

Corporate Office

Phone +1 610 225 6800, 550 E. Swedesford Road, Suite 400, Wayne, PA 19087, USA

Regional Offices

United States: Phone +1 919 544 8000, Toll Free 866 246 6990, cs@teleflex.com, 3015 Carrington Mill Boulevard, Morrisville, NC 27560, USA

Latin America: Phone +1 919 433 4999, la.cs@teleflex.com, 3015 Carrington Mill Boulevard, Morrisville, NC 27560, USA

International: Phone +353 (0)9 06 46 08 00, orders.intl@teleflex.com, Teleflex Medical Europe Ltd., IDA Business and Technology Park, Dublin Road, Athlone, Co Westmeath, Ireland

Teleflex, the Teleflex logo and Fortress are trademarks or registered trademarks of Teleflex Incorporated or its affiliates in the U.S. and/or other countries.

Information in this material is not a substitute for the product Instructions for Use. Not all products may be available in all countries. Please contact your local representative.

© 2025 Teleflex Incorporated. All rights reserved.

404774-EN · REV H · 10 25 PDF